



SIRIM Berhad
Industrial Biotechnology Research Centre, Building 19
Tel: 03-55446953/6960
Fax: 03-55446988

TEST REPORT

EVALUATION OF SKIN IRRITATION ON QUANTUM ION USING *IN VITRO* RECONSTRUCTED HUMAN EPIDERMAL MODEL EPIDERM™ SKIN IRRITATION TEST

Job No. J735/20

Report No. R735/20/B19/43

Sponsor:

Eva Energy Sdn Bhd,
12, Jalan Bandar 20,
Pusat Bandar Puchong,
47160 Puchong,
Selangor

Test Facility:

Industrial Biotechnology Research Centre (IBRC),
Building 19, SIRIM Berhad

Study Initiation Date:

29 June 2020

Experimental Start Date:

07 July 2020

Experimental End Date:

10 July 2020

Study Completion Date:

17 July 2020

TEST REPORT

REPORT NO:R735/20/B19/43

PAGE: 2 of 18

This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.

APPROVAL SIGNATURES

We, the undersigned, declare that the methods, results and data contained in this report faithfully reflect the procedures used and raw data collected throughout the study.



(SUZAINI BADRUDIN)

Reviewer

Industrial Biotechnology Research Centre

26 AUG 2020

Date



(NURHAYATI ARIFFIN)

Analyst

Industrial Biotechnology Research Centre

26 AUG 2020

Date



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 3 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

SUMMARY

EVALUATION OF SKIN IRRITATION ON QUANTUM ION USING *IN VITRO* EPIDERM™ SKIN IRRITATION TEST (SIT)

In vitro skin irritation test on Quantum Ion was performed according to the requirement of *In vitro* Skin Irritation: Reconstructed Human *Epidermis* Test Method –OECD Guidelines for Testing of Chemicals No. 439. This *in vitro* standard method was validated by European Centre of the Validation of Alternative Methods (ECVAM) as *in vitro* test method based on reconstructed human epidermis (RhE) technology. The test was conducted in line with Standard Operating Procedure (SOP) developed at MatTek Corporation. This test assesses irritability of both cosmetic ingredients and finished products.

The test was conducted to determine whether the test item cause irritation to the *in vitro* skin model EpiDerm™.

In vitro dermal irritation test consists of topical exposure of the Quantum Ion to reconstructed human epidermal model EpiDerm™ tissues, followed by a cell viability test. After 60 minutes of exposure, tissues were thoroughly rinsed, blotted to remove the test extract, and transferred to fresh medium. After a 24 hours incubation period, the medium was changed and tissues were incubated for another 18 hours. MTT [(3-4, 5 dimethyl triazole 2-yl) 2, 5-diphenyltetrazoliumbromide] assay was then performed by transferring the tissues to 6-well plates containing MTT medium (1 mg/mL). After 3 hours of incubation, the blue formazan salt formed by cellular mitochondria was extracted with 2.0 mL isopropanol / tissue. The optical density of the extracted formazan was determined using a spectrophotometer at 570 nm. Relative cell viability was calculated for each tissue as percentage (%) of the mean of the negative control tissues. The skin irritation potential was classified according to the remaining cell viability obtained after test item treatment.

The Quantum Ion did not reduce viability of the EpiDerm™ tissue to below 50 % of the negative control. Under the condition of this test, Quantum Ion is considered as **Non Irritant** to *in vitro* skin model EpiDerm™.



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 4 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

BACKGROUND

Skin irritation refers to reversible damage to the skin following the application of a test chemical for up to 4 hours as defined by the United Nations (UN) Globally Harmonized System of Classification and Labelling of Chemicals (GHS).

For chemical, the reconstructed human epidermal model was developed and designed to predict skin irritation potential of neat test substances in the context of identification and classification of skin irritation hazard according to the European Union (EU) classification system (R 38 or no label). Since the EU and GHS systems were harmonized in 2008, the procedure described in the SOP also allows for hazard identification of irritant substances in accordance to UN GHS.

In vitro skin irritation: Reconstructed Human Epidermis Test allows for assessment of irritation. A sufficient amount of extract was applied on the surface of the three dimensional reconstructed human epidermis (RhE). The RhE model is comprised of non-transformed human-derived epidermal keratinocytes, which have been cultured to form a multilayered, highly differentiated model of the human *epidermis*. It consists of organized basal, spinous and granular layers, and a multilayered *stratum corneum* containing intercellular lamellar lipid layers representing main lipid classes analogous to those found *in vivo*. After a certain incubation period, cell viability is assessed by MTT [(3-4, 5 dimethyl triazole 2-yl) 2, 5-diphenyltetrazoliumbromide] colorimetric test. The reduction of the viability of tissues exposed to chemicals in comparison to negative controls (treated with water) is used to predict the skin irritation potential.

Irritant chemicals are identified by their ability to decrease cell viability below defined threshold levels (*i.e.* $\leq 50\%$, for UN GHS Category 2). Depending on the regulatory framework and applicability of the Test Guideline, chemicals that produce cell viability above the defined threshold level may be considered non-irritants (*i.e.* $> 50\%$, No Category).



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 5 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

CONTENTS

APPROVAL SIGNATURES	2
SUMMARY	3
BACKGROUND	4
1.0 OBJECTIVE	6
2.0 STUDY TIMETABLE	6
3.0 MATERIALS	6
4.0 METHOD	9
5.0 DATA ANALYSIS	11
6.0 ACCEPTABILITY RANGES FOR QUALITY CONTROL	12
7.0 ACCEPTABILITY RANGES FOR NEGATIVE CONTROL	12
8.0 ACCEPTANCE CRITERIA FOR POSITIVE CONTROL	12
9.0 DATA INTERPRETATION PROCEDURE (PREDICTION MODEL)	13
10.0 RESULT AND DISCUSSION	13
11.0 CONCLUSION	15
12.0 RETENTION OF RECORDS AND TEST ITEM	15
13.0 REFERENCES	15
TABLES	
Table 1 Optical Densities (ODs)	16
Table 2 Blank Corrected Data	16
Table 3 Classification of Test Item	17
FIGURES	
Figure 1 Relative Viability	18



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 6 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

1.0 OBJECTIVE

The objective of this study is to predict skin irritation potential of test item using *in vitro* skin model EpiDerm™.

2.0 STUDY TIMETABLE

- 2.1 **Receipt of reconstructed human epidermal model EpiDerm™ tissues:**
07 July 2020
- 2.2 **Tissue Conditioning**
07 July 2020
- 2.3 **Pre-Incubation**
08 July 2020
- 2.4 **Treatment**
08 July 2020
- 2.5 **Change Medium**
09 July 2020
- 2.6 **MTT Viability Test**
10 July 2020
- 2.7 **Optical Density Reading**
10 July 2020
- 2.8 **Data Analysis**
10 July 2020 – 17 July 2020

3.0 MATERIALS

3.1 Test Item

- 3.1.1 Test item: Quantum Ion
- 3.1.2 Sample marking: Copper Ion
- 3.1.3 Date received: 29 June 2020
- 3.1.4 Physical appearance: Liquid
- 3.1.5 Colour: Clear blue
- 3.1.6 Physical Chemical Properties Data: Not provided
- 3.1.7 Quantity received : 250 ml + 50 ml + 50 ml
- 3.1.8 pH: Not provided



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 7 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

3.1.9 Storage condition: Room temperature

3.1.10 Solubility: Not provided

3.1.11 Stability: Not provided

3.1.12 Expiration date: Not provided

3.2 Test System

3.2.1 Test System: Reconstructed human epidermal model EpiDerm™

The reconstructed human epidermal model EpiDerm™ (EPI-200, MatTek, Ashland, USA) consists of normal human-derived epidermal keratinocytes, which have been cultured to form a multilayered highly differentiated model of the human epidermis. It consists of organized basal, spinous and granular layers, and a multilayered stratum corneum containing intercellular lamellar lipid layers arranged in patterns analogous to those found *in vivo*.

The EpiDerm™ tissues are cultured on specially prepared cell culture inserts, containing 24 tissues using serum free medium. Ultrastructurally, the EpiDerm Skin Model closely parallels human skin, thus providing a useful *in vitro* means to assess dermal irritancy and toxicology.

3.2.1.1 Lot No: 33066

3.2.1.2 Production Date: 02 July 2020

3.2.1.3 Date of Shipping: 03 July 2020

3.2.1.4 Receipt of EpiDerm™:07 July 2020

3.2.1.5 Visual quality control of the skin: All tissues in good condition

3.2.1.6 The EpiDerm™ System is manufactured according to defined quality assurance procedures. All biological components of the epidermis and the culture medium are tested by manufacturer for viral, bacterial, fungal and mycoplasma contamination. MatTek determines the ET-50 value following exposure to Triton X-100 (1%) for each EpiDerm™ lot. The ET-50 must fall within a range established based on a historical database of results or acceptability ranges for quality control based on OECD Test Guidelines. The quality control value is presented in 10.3.

3.3 Reagent

3.3.1 Assay Medium: EPI-100-NMM-SIT / Assay Medium

3.3.1.1 Lot No. : 121119TVKD

3.3.1.2 Sterility : Sterile

3.3.1.3 Expiration Date : 01/08/2020

3.3.1.4 Storage : Refrigerator (5 ± 3 °C)

3.3.1.5 Manufacturer : Mattek Corporation



MS ISO/IEC 17025
TESTING
Samm No.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 8 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

- 3.3.2 Phosphate Buffered Saline without Calcium and Magnesium
Lot No.: 111219RAHA
Expiration Date : 12/12/2020
Manufacturer : Mattek Corporation
- 3.3.3 Phosphate Buffered Saline without Calcium and Magnesium, contain the following:
- 3.3.3.1 Sodium Chloride 8 g/L
Lot No.: K50722004912
Manufacturer : MERCK
- 3.3.3.2 Potassium chloride 0.2 g/L
Lot No.: 1409BI4J35501
Manufacturer : Bio Basic Canada Inc
- 3.3.3.3 Anhydrous potassium dihydrogen orthophosphate 0.2 g/L
Lot No.: 1211ACN1K12062801
Manufacturer : Bio Basic Canada Inc
- 3.3.3.4 Anhydrous disodium hydrogen orthophosphate 1.15 g/L
Lot No.: 1306ACK2NA12020101
Manufacturer : Bio Basic Canada Inc
- 3.3.4 MTT – 2mL (5mg/mL)
Lot No. : MKBW0025V
CAS No. : 298-93-1
Manufacturer : SIGMA
- 3.3.5 MTT Diluent – 8mL
As above 3.3.1
- 3.3.6 Extractant Solution – Isopropanol
Lot No. : 632261
Storage : Room Temperature
CAS No. : 67-63-0
Manufacturer : Fisher Scientific
- 3.3.7 Positive Control: 5 % SDS Solution
Part No. : TC-SDS-5 %
Lot No. : 121219BBA
Expiration Date : 12/12/2020
Storage : Room Temperature
Manufacturer : Mattek Corporation
- 3.3.8 Negative Control: As above 3.3.2

3.4 Material

- 3.4.1 Nylon Mesh (EPI-MESH)
3.4.2 Lot no. : 0551428-00
3.4.3 Expiry Date : 31/12/2029
3.4.4 Manufacturer : MatTek Corporation



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 9 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

4.0 METHOD

Upon receipt of the reconstructed human epidermal model EpiDerm™, the tissue kits and the solutions were stored according to the manufacturer's directions for unpacking and storage.

4.1 Tests for Interference of Test Item

Tissue-Binding of Coloured or Staining Materials

30 µL of test item was applied onto three single EpiDerm™ tissues. In parallel, a tissue was exposed to negative control. After 60 minutes of exposure, tissues were thoroughly rinsed, blotted to remove the test extract, and transferred to fresh medium.

After a 24 hours incubation period, the medium was changed and tissues were incubated for a further 18 hours. The medium was then changed again and incubated for 3 hours. After incubation was completed, the tissues were rinsed and extracted using 2.0 mL of isopropanol per tissue. The optical density of the extracted tissue was determined using spectrophotometer at 570 nm.

4.1.2 Test for Interference of Test Item with MTT

30 µL of test item of the test item was added to 1 mL of the 1 mg/mL MTT and incubated at (37 ± 1)° C, (5 ± 1)% CO₂, 95 % relative humidity for 60 minutes. Untreated MTT medium was used as control. If the MTT solution turns blue/purple, the test item reduces MTT and additional functional check must be performed.

4.1.3 Test for Mesh Compatibility

A mesh was placed on a glass slide and 30 µL of test item was applied. After 60 minutes of exposure, the reaction between the mesh and test item was observed under microscope.

4.2 Tissue Conditioning – Day 0

Under sterile conditions, a sealed 24-well plate containing the EpiDerm™, was opened. Visual inspection of each insert containing the epidermal tissue was done prior to tissue conditioning.

0.9 mL of assay medium was dispensed into each well of six-well plate, the EpiDerm™ tissue cultures were transferred into the wells and the plates was incubated for (60 ± 5) minutes at (37 ± 1)° C, (5 ± 1)% CO₂, 95 % relative humidity. At the end of the first 60 minutes pre-incubation, the assay medium was renewed for further pre-incubation. Each insert was aseptically transferred into well containing 0.9 mL assay medium and pre-incubation was done at (37 ± 1)° C, (5 ± 1)% CO₂, 95 % relative humidity for (18 ± 3) hours.



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 10 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

4.3 Dosing protocol – Day 1

30 μL of test item was applied onto three single EpiDerm™ tissues. Negative and positive controls were conducted in parallel using identical method to the dosed cultures. The cultures were incubated at $(37 \pm 1)^\circ\text{C}$, $(5 \pm 1)\%$ CO₂, 95 % relative humidity for (35 ± 1) minute.

After incubation, each dosed EpiDerm™ tissue was removed from the incubator and placed at room temperature in the biological safety cabinet for 25 minutes. At the end of the exposure period, the EpiDerm™ tissue cultures were removed from the assay plates and gently rinsed with phosphate buffered saline to eliminate any residual test material. The EpiDerm™ tissue cultures were then transferred into wells containing 0.9 mL assay medium and the plates incubated at $(37 \pm 1)^\circ\text{C}$, $(5 \pm 1)\%$ CO₂, 95 % relative humidity for (24 ± 2) hours.

4.4 Change Medium – Day 2

After 24 hour's incubation, the medium was changed and EpiDerm™ tissues were incubated for a further 18 hours.

4.5 MTT Viability Test

4.5.1 Day 3

The MTT assay was performed by transferring each EpiDerm™ tissues to 6-well plates containing 300 μL MTT medium (1 mg/mL) in each well. After 3 hours of MTT incubation, the blue formazan salt formed by cellular mitochondria was extracted with 2.0 mL isopropanol / tissue (extractant solution). The extraction plates were sealed with parafilm and agitated for at least 2 hours at room temperature. As alternative, overnight extraction at room temperature in the dark also possible.

4.5.2 At the end of the extraction period, EpiDerm™ tissue was pierced with an injection needle and the extract was decanted into the well from which insert was taken. The insert was then discarded. The extraction solution was pipetted up and down to ensure complete mixing. Finally, 200 μL was transferred into a 96 well microtiter plate for absorbance measurement (OD=optical density) at 570 nm without using a reference filter. 200 μL of isopropanol was used as blank.

Relative cell viability for each tissue was calculated as percentage (%) relative to mean of the negative control tissues viability.



TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 11 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

5.0 DATA ANALYSIS

5.1 Tests for Interference of Test Item

5.1.1 Assessment of Coloured or Staining Materials

The percentage difference of the optical density between a coloured test item (CTI) and the negative control (NC) was calculated according to the following formula:

$$[(\text{Mean OD}_{\text{CTI}} - \text{Mean OD}_{\text{NC}}) / \text{Mean of OD}_{\text{NC}}] \times 100$$

OD reading of treated tissue by coloured test item	Action
Below 5 % (<5 %) of the negative control ; Tissue viability determined in MTT Assay not close to classification cut-off (50%)	Correction of the results is not necessary
Between 5 % and 30 % of the negative control	Further test on more tissue
Above 30 % (>30 %) of the negative control	Additional step and expert judgment ; Incompatible with the test system

The real MTT OD (unaffected by interference with the color or staining materials) was calculated using the following formula:

$$\text{OD} = \text{OD Coloured tissue (MTT assay)} - \text{OD Coloured tissue (no MTT assay)}$$

5.1.2 Test for Interference of Test Item with MTT

The test item is presumed to have reduced the MTT if the MTT solution colour turns blue/ purple.

5.1.3 Test for Mesh Compatibility

If there is reaction between test item and the mesh, the test item has to be applied without using mesh as a spreading aid.

5.2 Workbook EpiDerm™-SIT

A blank, password protected MS EXCEL workbook EpiDerm™-SIT-SPREAD.XLS was provided by MatTek Corporation. The workbook consists of two single spreadsheets named: Import and Spread.

5.3 Raw Data of Optical Densities (ODs)

Raw data of optical densities generated by the microplate reader (without blank subtraction) were copied from the reader software and then pasted into the Import spreadsheet of the Excel workbook. The blank corrections, calculation of results and statistical parameters are automatically calculated in the Spread spreadsheet of the workbook.



TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 12 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

5.4 Calculation

- 5.4.1 After data entry, the spreadsheet performs the following calculations:
- 5.4.2 Blank correction
- 5.4.3 For each individual tissue treated with a test item (TI), the positive control (PC) and the negative control (NC) the individual relative tissue viability was calculated according to the following formulas;

$$\text{Relative viability TI (\%)} = [\text{OD}_{\text{TI}} / \text{Mean of OD}_{\text{NC}}] \times 100$$

$$\text{Relative viability NC (\%)} = [\text{OD}_{\text{NC}} / \text{mean of OD}_{\text{NC}}] \times 100$$

$$\text{Relative viability PC (\%)} = [\text{OD}_{\text{PC}} / \text{mean of OD}_{\text{NC}}] \times 100$$

- 5.4.4 For each test item, negative control and positive control, the mean relative viability of the three individual tissues was calculated and used for classification according to the Prediction Model (Refer to 8.0).

6.0 ACCEPTABILITY RANGES FOR QUALITY CONTROL

	Lower acceptance limit	Upper acceptance limit
EpiDerm™ SIT (EPI-200) (1% Triton X-100)	$\text{ET}_{50} = 4.0 \text{ hour}$	$\text{ET}_{50} = 8.7 \text{ hour}$

7.0 ACCEPTABILITY RANGES FOR NEGATIVE CONTROL OD VALUES OF THE TEST METHODS

	Lower acceptance limit	Upper acceptance limit
EpiDerm™ SIT (EPI -200)	≥ 1.0	≤ 2.5

8.0 ACCEPTANCE CRITERIA FOR POSITIVE CONTROL

The assay meets the acceptance criterion if the mean viability of positive control tissues expressed as percent of the negative control tissues is $\leq 20\%$. The standard deviation shall be below 18 for all substances and controls.

	Mean of Viability (%)	Standard Deviation of Viability (SD %)
EpiDerm™ SIT (EPI-200)	≤ 20	<18



TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 13 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

9.0 DATA INTERPRETATION PROCEDURE (PREDICTION MODEL)

An irritant is predicted if the mean relative tissue viability of three individual tissues exposed to the test substance is reduced below 50% of the mean viability of the negative controls.

<i>In vitro</i> result	<i>In vivo</i> prediction
mean tissue viability ≤ 50%	Irritant (I), (R38 or GHS category 2)
mean tissue viability > 50%	non-irritant (NI)

10.0 RESULT AND DISCUSSION

10.1 Tests for Interference of Test Item

10.1.1 Assessment of Non-Coloured or Staining Materials

There was no colour change in deionized water, therefore the test item did not have potential to stain the tissue.

10.1.2 Test for Interference of Test Item with MTT

There was no change in MTT colour therefore the test item did interfere with MTT.

10.1.3 Test for Mesh Compatibility

There was no reaction between the mesh and the test item; therefore the test item is compatible to the mesh.

10.2 Viability Measurement

10.2.1 Raw data of optical densities (ODs) of Blank

Refer To Table 1

10.2.2 Blank Correction

Refer to Table 2

10.3 The Quality Control value meets the acceptance range criteria

EpiDerm™ SIT (EPI-200)	Hour
ET ₅₀	6.12

10.4 The negative control OD value meets the acceptance range criteria

EpiDerm™ SIT (EPI-200)	Lower OD	Upper OD
Negative Control	1.983	2.208



MS ISO/IEC 17025
TESTING
SAMM NO.676



TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 14 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

10.5 The positive control OD value meets the acceptance range criteria

EpiDerm™ SIT (EPI-200)	Mean of Viability (%)	Standard Deviation of Viability (SD %)
Positive Control	7.3	0.16

10.6 Classification

Test item, negative control and positive control are qualified according to prediction model. The mean relative viability of the three individual tissues was calculated and used for classification

Refer to Table 3

10.7 Graph

10.7.1 The spreadsheet shows a graph of the results (% of relative viability ± standard deviation)

Refer to Figure 1

10.8 Prediction Model

	Mean of Viability (%)	Standard Deviation of Viability (SD %)	In vitro result	In vivo prediction
Quantum Ion	93.8	4.69	mean tissue viability > 50 %	Non-Irritant (NI)
Negative Control	100.0	5.39	mean tissue viability > 50 %	Non-Irritant (NI)
Positive Control	7.3	0.16	mean tissue viability ≤ 50 %	Irritant (I), (R38 or GHS category 2)



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 15 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

11.0 CONCLUSION

Under the condition of this test, Quantum Ion is considered as **Non-Irritant** *in vitro* skin model EpiDerm™.

12.0 RETENTION OF RECORDS AND TEST ITEM

One report will be forwarded to the Sponsor. The other report, together with all generated raw data is maintained at the Industrial Biotechnology Research Centre Archives.

13.0 REFERENCES

- 13.1 OECD (2013), In Vitro Skin Irritation: Reconstructed Human *Epidermis* Test Method – OECD Guidelines for Testing of Chemicals No. 439
- 13.2 Protocol for: In Vitro EpiDerm™ Skin Irritation Test (EPI-200-SIT) – Reconstructed Human Epidermal Model EpiDerm (EPI-200-SIT). For use with MatTek Corporation
- 13.3 (LWI-238-43) : In Vitro EpiDerm™ Skin Irritation Test (EPI-200-SIT)



MS ISO/IEC 17025
TESTING
SAMM NO.676

TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 16 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

Table 1 Optical Densities (ODs) of Blank

Optical Densities (ODs) of Blank	Mean Optical Densities (ODs) of Blank
0.0387	
0.0459	
0.0418	
0.0414	
0.0400	
0.0403	0.0414

Table 2 Blank Corrected Data

	Tissue	Raw data		Blank corrected data		Mean	% of Viability
		1	2	1	2		
Quantum Ion	1	1.8736	1.9095	1.832	1.868	1.850	88.5
	2	2.0746	2.0841	2.033	2.043	2.038	97.5
	3	2.0634	2.0047	2.022	1.963	1.993	95.3
Negative Control	1	2.2734	2.2253	2.232	2.184	2.208	105.6
	2	2.1604	2.0834	2.119	2.042	2.081	99.5
	3	1.9548	2.0948	1.913	2.053	1.983	94.9
Positive Control	1	0.1983	0.1975	0.157	0.156	0.157	7.5
	2	0.1916	0.1914	0.150	0.150	0.150	7.2
	3	0.1955	0.1908	0.154	0.149	0.152	7.3



TEST REPORT

REPORT NO:R735/20/B19/43	PAGE: 17 of 18
This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.	

Table 3 Classification of Test Item

	Mean of OD	SD of OD	Mean of Viability [%]	SD of Viability	CV [%]	<i>In vitro</i> result	Classification
Quantum Ion	1.960	0.098	93.8	4.69	5.00	Non-Irritant	Qualified
Negative Control	2.091	0.113	100.0	5.39	5.39	Non-Irritant	Qualified
Positive Control	0.153	0.003	7.3	0.16	2.17	Irritant	Qualified

OD- Optical Density, SD- Standard Deviation, CV- Coefficient of Variation



TEST REPORT

REPORT NO:R735/20/B19/43

PAGE: 18 of 18

This report is NOT a Quality Assurance Certificate NOR an Approval Permit. This report refers only to samples submitted by the customer to SIRIM Berhad and tested by SIRIM Berhad. This report shall not be reproduced, except in full and shall not be used for advertising purposes by any means or forms without written approval from President & Chief Executive of SIRIM Berhad.

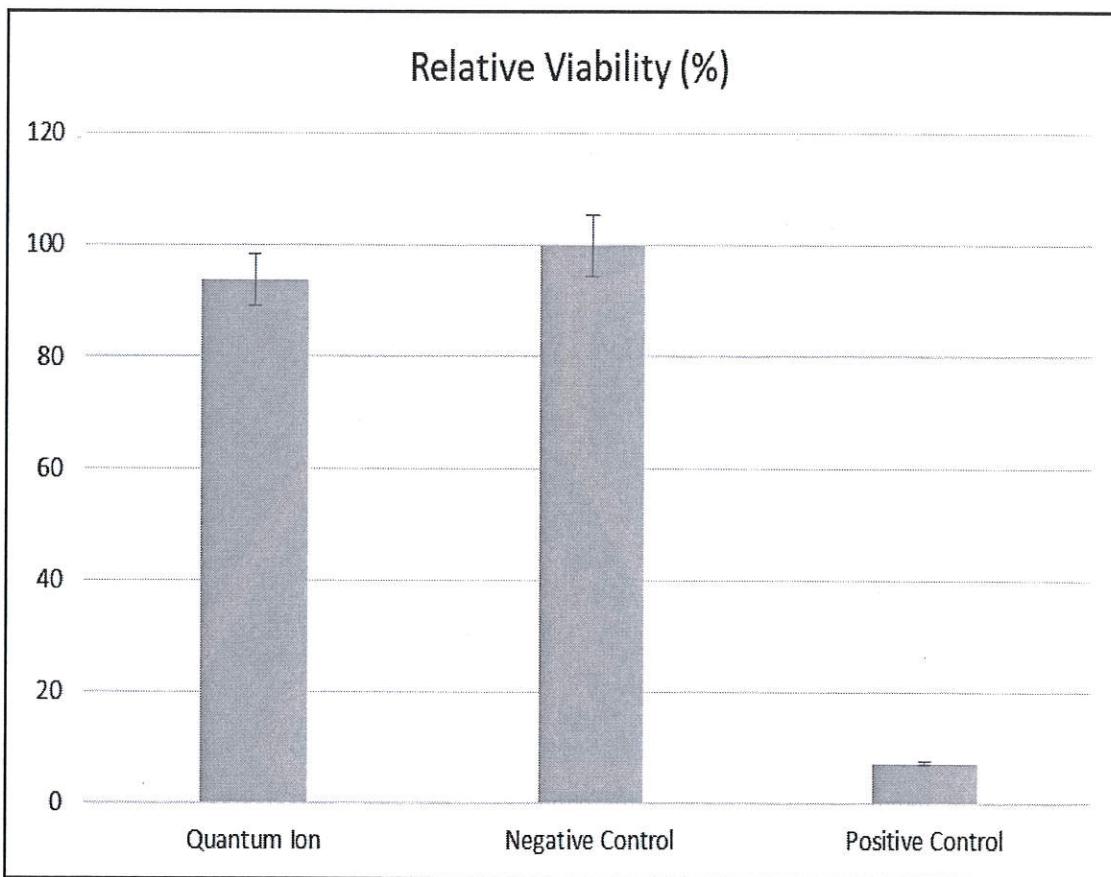


Figure 1 Relative Viability



MS ISO/IEC 17025
TESTING
SAMM NO.676

CONDITIONS RELATING TO THE USE OF SIRIM BERHAD TEST REPORT

1. A Test Report will be issued in respect of Testing Services conducted and shall relate only to the sample actually tested. SIRIM Berhad makes no warranty whatsoever and the Applicant shall not represent in any manner that any duplication or mass production of the Product is same as the Sample actually tested or that SIRIM Berhad has tested any of the duplicated or mass produced Product.
2. The Test Report shall not be amended, changed, varied or modified in any manner whatsoever by the Applicant or otherwise.
3. If the Test Report is to be furnished to any third party or to the public, each such Test Report shall be furnished in full, legible and in its entirety.
4. The Test Report shall not be reproduced and shall not in any event be used for any advertising purposes or whatsoever without written approval from the President & Chief Executive of SIRIM Berhad of No. 1, Persiaran Dato' Menteri, Building 5, Section 2, P. O. Box 7035, 40700 Shah Alam, Selangor Darul Ehsan.
5. Customer (Applicant/Manufacturer/Factory, etc.) is not permitted to use any SIRIM Berhad, other SIRIM Berhad's subsidiaries logo on packaging, sample's manual, technical specification, brochures/flyers or any other means.
6. If such approval is obtained from the President & Chief Executive, the Applicant may only include the phrase, "A sample of this product has been tested by SIRIM Berhad ... (Test Report No) ... (dated) ... (for what test) ... (to which standard)" or such similar words which stress that only the Sample was actually tested. This phrase shall only be used for the purpose of product advertisement or product promotion (eg; brochures). For avoidance of doubt, the statement shall not be used on the sample and packaging of the sample.
7. In the event there is an investigation from a Government Regulatory Agency concerning the applicant's Test Report, SIRIM Berhad may disclose the information pertaining to the Test Report for purposes of such investigation.
8. Further or in the alternative, it is strictly forbidden to represent in any manner whatsoever that SIRIM Berhad and/or other SIRIM's subsidiaries has endorsed, approved or validated the Product of the Applicant in any manner whatsoever.
9. In the event the applicant is found in breach of this provision, SIRIM Berhad and/or other SIRIM's subsidiaries without prejudice to any other rights and remedies may take whatever action necessary including but not limited to:
 - a) Informing and placing a notice in the media;
 - b) Obtaining an injunction from Court (cost on a solicitor-client basis to be borne by the Applicant);
 - c) Refusing to accept any further Product for Testing Services from the Applicant or whatsoever related to the Applicant, whether subsidiary or otherwise;
 - d) Instructing the Applicant to withdraw and recall the advertisement, statement or document in question and advertise a clarification and apology to SIRIM Berhad and/or other SIRIM's subsidiaries twice in a national publication of SIRIM Berhad's choice at the Applicant's sole cost; and
 - e) Informing or lodging a report pertaining the Applicant's Test Report with the relevant authorities.
10. Certified true copies of the Test Report may be issued upon request by the applicant upon payment of the relevant fee.
11. Corrections to test report shall only be allowed within 6 months from issuance date of the Test Report of the relevant fee and shall be limited to maximum 3 times, after either case whichever occurs earlier, a new Test Report shall be issued and replace the previous one (having error(s) or lack of information) with relevant fee. Issuance of Supplementary Report to the original Test Report shall be for the followings:
 - a) Misprints and typo errors;
 - b) Missing technical information;
 - c) Test data not reported;
 - d) Mistake in reporting of test data.
12. Any amendment requested from customers on the test report issued shall be in writing.
13. SIRIM reserves the right in its sole discretion to terminate or modify this permission.